Oral-B®

Professional Care WaterJet

MD 15 A

powered
by **BRAUN**

Type 3718











Deutsch	4, 36, 41
English	6, 36, 41
Français	8, 36, 41
Español	10, 36, 41
Português	12, 37, 41
Italiano	14, 37, 41
Nederlands	16, 37, 41
Dansk	18, 38, 41
Norsk	20, 38, 41
Svenska	22, 38, 41
Suomi	24, 38, 41
Ελληνικά	26, 39, 41
Türkçe	29, 41
Русский	31, 39, 41
عربي	35, 33, 41

Internet:

www.oralb.com
www.braun.com

- (D)** Haben Sie Fragen zu diesem Produkt?
(A) Rufen Sie an (gebührenfrei):
00 800 27 28 64 63
- (CH)** **08 44 - 88 40 10**
- (GB)** If you should require further assistance, please call our Consumer Service line.
0 800 73 11 792 (free call)
- (IRL)** **Helpline**
1 800 509 448 (free call)
- (F)** Appelez le Service Consommateurs Oral-B :
0 810 365 855 (numéro azur)
- (B)** **02-711 92 11**
- (LU)** Vous avez des questions sur ce produit?
Appelez:
4 00 50 51
- (E)** Servicio al consumidor para España:
901 11 61 84
- (P)** Serviço ao Consumidor para Portugal:
808 20 00 33
- (I)** Servizio consumatori:
02-667 86 23
- (NL)** Heeft u vragen over dit product?
Bel Braun/Oral-B Consumenten-lijn:
0800-4 45 53 88
- (DK)** Har du spørgsmål om produktet? Så ring
70 15 00 13
- (N)** Spørsmål om dette produktet? Ring
22-72-88-10
- (S)** Frågor om apparaten? Ring Kundservice
020 - 21 33 21
- (FIN)** Onko Sinulla kysyttävää tuotteesta? Soita
09 452 872
- (GR)** **Τηλέφωνο για πληροφορίες**
210 9478700
- (TR)** **Danışma Hattı**
0212-473-75-85
- (RUS)** **Телефонный помощь**
8-800-200-1444

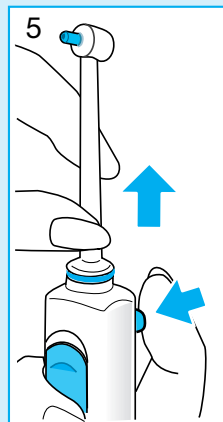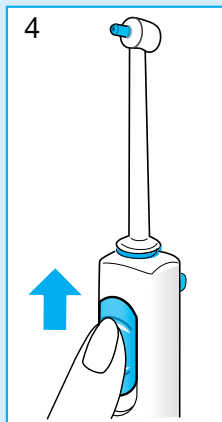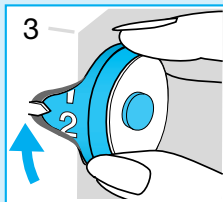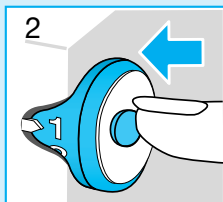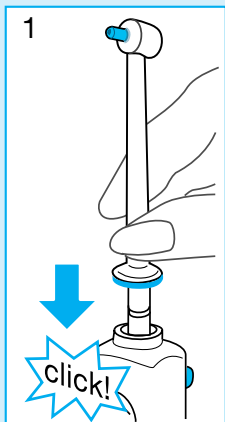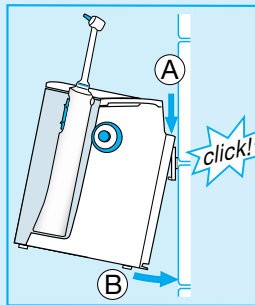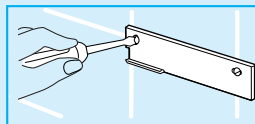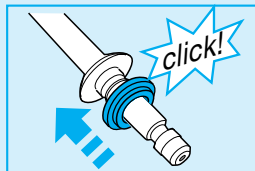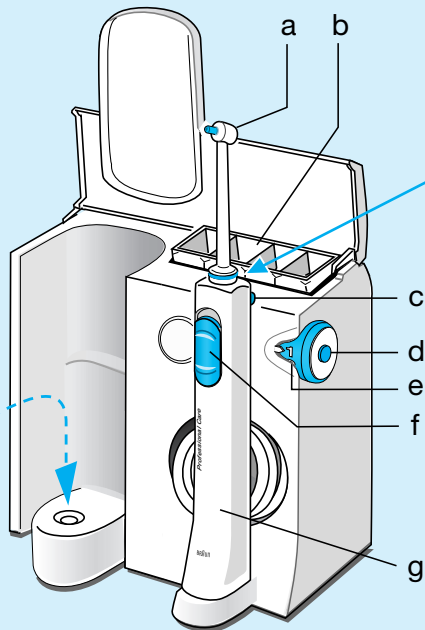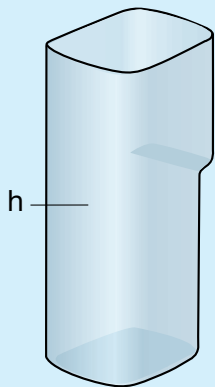Deutsch

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English

The Oral-B Professional Care WaterJet removes remaining food particles from between your teeth, while it also effectively massages your gums. Always use the WaterJet after brushing your teeth.

Please read the use instructions carefully before first use.

Important

The appliance must be mounted in a way that it cannot fall into water.

The plug should always be accessible. Periodically check the cord for damage. A damaged or non-functioning appliance should no longer be used. Take it to a Braun Service Centre.

Description

- Ⓐ Irrigator jet
- Ⓑ Jet compartment
- Ⓒ Jet release button
- Ⓓ On/off switch
- Ⓔ Water pressure control dial
- Ⓕ Handle switch (water stop and go)
- Ⓖ Irrigator handle
- Ⓗ Water container (600 ml/20 fl. oz.)

Specifications

Power supply: see bottom of housing

Power intake: 24 W

Using the Irrigator


Before using for the first time, fill the water container and allow the water to run through completely by holding the irrigator handle over a washbasin.

1. Attach one of the irrigator jets Ⓐ to the handle Ⓖ until it locks in place. Fill the water container with water and add mouthwash, if desired.
2. Turn on the irrigator by pressing the on/off switch Ⓓ. (When doing so, the handle switch Ⓕ should be in the lower position = water stop.)
3. Select the water pressure with the pressure control dial Ⓔ (1 = gentle, 5 = strong). We recommend that you start with setting (1), especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over the washbasin and guide the jet to your teeth and gums. Slide the handle switch up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.

The irrigator jet offers two different functions depending on the position of the jet tip:

Single jet:

Jet tip in position (a) to clean the teeth and the interdental spaces and remove resistant food particles.

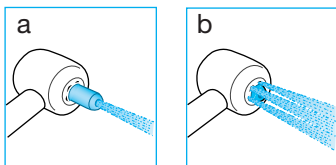Multi jet:

Jet tip in position (b) to massage the gums and keep them healthy.

First, turn off the water stream by sliding the handle switch down (to avoid unintended water spraying).

The jet can easily be switched from one function to the other by pressing the jet tip in or out.

Begin with the single jet. Guide the water stream to the teeth and especially to the interdental spaces.

Then massage the gums with the multi jet.

To ensure the jet functions, a small sieve is incorporated in the lower end of the shaft. Should you notice a reduced water pressure, check the sieve for any deposits and remove them, e.g. by carefully tapping the shaft on a flat surface.

Max. operating time: 15 minutes

Cooling off time: 2 hours

After use

Always empty the water container entirely to prevent bacteria from collecting in the standing water. If your dentist has recommended

use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use to prevent any deposit build-up.

Press the on/off switch to turn off the irrigator.

To take off the irrigator jet, press the jet release button. The jets can be stored in the jet compartment. Wipe the handle dry and place it back onto its base.

Cleaning

Regularly clean the housing with a damp cloth.

The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

This product conforms to the European Directives EMC 89/336/EEC and Low Voltage 73/23/EEC.

At the end of the product's useful life, please dispose of it at appropriate collection points provided in your country.

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English

Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send

the complete appliance with your sales receipt to an authorised Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

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Hong Kong

Audio Supplies Company,
Room 506,
St. George Building,
2 ICE House Street,
Hong Kong,
☎ (852) 2524 9377

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Tehran Bouran Company
145 Ghaem Magham Farahani
Tehran
☎ +982 18 31 27 66

Ireland (Republic of)

Gavin's Electronics,
83-84, Lower Camden Street,
Dublin 2,
☎ 1800 509 448

Israel

S. Schestowitz Ltd.,
8 Shacham Str.,
Tel-Aviv, 49517,
☎ 1 800 335 959

Italia

Servizio Consumatori Braun
Gillette Group Italy S.p.A.,
Via G. B. Pirelli, 18,
20124 Milano,
☎ 02 / 6678623

Jordan

Interbrands
Wadi Al-Sir, Adwiat alhekma
St. Mekadabi bld.
Amman
☎ +692 582 75 67

Kenya

Radbone-Clark Kenya Ltd.,
P.O. Box 40833,
Mombasa Road,
Nairobi,
☎ 2 82 12 76

Korea

Gillette Korea Ltd.
144-27 Samsung-dong,
Kangnam-ku,
Seoul, Korea,
☎ 080-920-6000

Kuwait

Union Trading Company,
Braun Service Center,
P. O. Box 28 Safat,
Safat Code 13001, Kuwait,
☎ 04 83 32 74

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Philippines

Gillette Philippines Inc.,
Corporate Corner Commerce Avenues
20/F Tower 1, IL Corporate Centre
1770 Muntinlupa city
☎ 02-771071 02-06/-16

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South Africa (Republic of)

Fixnet After Sales Service,
159 Queen Street, Kensington South,
P.O. Box 751770,
Johannesburg 2094,
Kensington South,
☎ 11 615 6765

St. Maarten

Rupchand Sons n.v. (ram's),
P.O. Box 79,
Philipsburg,
Netherlands Antilles,
☎ 52 29 31

St. Thomas

Boolchand's Ltd.,
31 Main Street,
P.O. Box 5667,
00803 St. Thomas,
US Virgin Islands,
☎ 340 776 0302

Suomi

Gillette Group Finland Oy,
P.O. Box 9,
Niittykatu 8, PL 9,
02200 Espoo,
☎ 09-45 28 71

Sverige

Gillette Group Sverige AB,
Dept. Sweden,
Stockholm Gillette
Räsundavägen 12,
Box 702,
16927 Solna,
☎ 020-21 33 21

Syria

Ahmed Hadaya Company
Hadaya building
Ain Keresh
Unisyria, P.O. Box 35002,
Damascus,
☎ 963 011-231433

Taiwan

Audio & Electr. Supplies Ltd.,
Brothers Bldg., 10th Floor,
85 Chung Shan N Rd., Sec. 1,
Taipei (104),
☎ (886) 02 2523 3283

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Yemen (Republic of)

Saba Stores for Trading,
26th September Street,
P.O. Box 5278,
Taiz,
☎ 4-25 23 88

Yugoslavia

BG Elektronik,
Bulevar kralja Aleksandra 34,
11000 Beograd,
☎ 11 3240 030

Oral-B Professional Care WaterJet MD 15 A

Garantiekarte
Guarantee Card
Carte de garantie
Tarjeta de garantía
Cartão de garantia
Carta di garanzia
Garantiebewijs
Garantibevis
Köpbevis
Takuukortti
Κάρτα εγγυήσεως
Гарантийный талон

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